CHLAMYDIA CEL PN

English
Product Code: KC3

FIGURE 1: CHLAMYDIA CEL PN DIAGRAM FOR USE

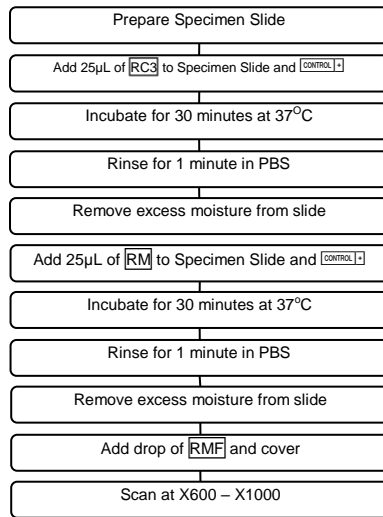
TABLE 1: SENSITIVITY, SPECIFICITY, & OTHER DATA ON THE CHLAMYDIA CEL PN
 TABLEAU 1: SENSIBILITÉ, SPÉCIFICITÉ ET AUTRES DONNÉES DU TEST CHLAMYDIA CEL PN
 TABELLE 1: SENSIBILITÄT, SPEZIFIZITÄT UND ANDERE DATEN ZUM CHLAMYDIA CEL PN
 TABELLA 1: SENSIBILITA', SPECIFICITA' ED ALTRI DATI SULLA CHLAMYDIA CEL PN
 TABLA 1: SENSIBILIDAD, ESPECIFICIDAD Y OTROS DATOS DEL CHLAMYDIA CEL PN
 TABELA 1: SENSIBILIDADE, ESPECIFICIDADE, & OUTROS DADOS DO CHLAMYDIA CEL PN
 TABELA 1: CZUŁOŚĆ, SWOISTOŚĆ I INNE DANE NA TEMAT ZESTAWU CHLAMYDIA CEL PN

Trial Essai Versuch Prova Prueba Teste Próba	Sensitivity Sensibilité Sensitivität Sensibilita' Sensibilidad Sensibilidade Czułość	Specificity Spécificité Spezifität Specificita' Especificidad Especificidade Swoistość	Repeatability Répétabilité Wiederholpräzision Repetibilita' Repetibilidad Repetição Powtarzalność	Reproducibility Reproductibilité Reproduzierbarkeit Riproducibilita' Reproducibilidad Reprodutibilidade Odtwarzalność
A	100%	100%	-	-
B	-	-	100% Correlation	100% Correlation
Not cross reactive with / Pas de Réaction Croisée avec / Keine Kreuzreaktionen mit / Non mostra reazione crociata con / No muestra reacción cruzada con / Sem reacção cruzada com / Brak reakcji krzyżowej z: <i>Staphylococcus epidermidis, Staphylococcus aureus, Beta haemolytic Streptococcus group A, Streptococcus pneumoniae, Haemophilus influenzae, Haemophilus parainfluenzae, Escherichia coli, Klebsiella pneumoniae, Proteus mirabilis, Pseudomonas aeruginosa, Bacteroides fragilis, Viridans Streptococcus, Candida albicans, Mycoplasma pneumoniae, Chlamydia trachomatis, Chlamydia psittaci, Respiratory Syncytial Virus, Adenovirus, Parainfluenzae type 3, Influenzae type B, Cytomegalovirus</i>				

EXPLANATION OF SYMBOLS

- Consult Instructions for Use
- In Vitro Diagnostic Medical Device
- Temperature Limitation
8°C
2°C
- Batch
- Control Positive
- Use By/Expiration Date
- Do Not Re-use

Cellabs Pty Ltd
 Unit 7, 27 Dale Street
 Brookvale, NSW 2100 Australia
 Tel: +61 2 9905 0133 Fax: +61 2 9905 6426
 Web: <http://www.cellabs.com.au>
 Email: sales@cellabs.com.au

WMDE
 Bergerweg 18
 6085 AT Horn
 The Netherlands

en fr de it es pt pl
 LC3.15.1
 23 May 2017

CE 0843

INTENDED USE AND PRINCIPLE OF THE TEST

The Chlamydia Cel Pn If Test is an *in vitro* indirect immunofluorescence test for the direct detection of *Chlamydia pneumoniae* organisms in clinical specimens. A monoclonal antibody to *C. pneumoniae* binds specifically to the organisms present in fixed specimens and a second antibody, FITC-conjugated goat anti-mouse Ig, stains the *C. pneumoniae* organisms.

CONTENTS OF THE KIT

		KC3 Standard	Bulk
	Chlamydia Cel Pn Mab Reagent	1.25mL	5mL
	Anti-mouse Ig-FITC Reagent	1.25mL	5mL
	Positive Control Slide (Single Use only)	1	-
	Mounting Fluid	2.5mL	-
	Tests	50	200

Materials are supplied ready for use. Store at 2-8°C. Expiry dates are clearly marked on each kit component and the box. Expiry dates do not change once opened.

MATERIALS REQUIRED BUT NOT PROVIDED

Microscope slides with 6-8mm diameter wells; precision pipette for delivering 25µL; acetone or methanol for specimen fixation; humid chamber; wash bath; phosphate buffered saline (PBS) for washing step; cover slips; non-fluorescing immersion oil; and fluorescence microscope with filter system for FITC (maximum excitation wavelength 490nm, mean emission wavelength 530nm) and x600-x1000 magnification.

PRECAUTIONS

For *in vitro* diagnostic use only. Do not use after the expiry date shown on the label. If protective packaging is damaged, contact your local distributor and ask for a replacement. Do not mix components from different kits. The Chlamydia Cel Pn Mab Reagent and Anti-mouse Ig-FITC Reagent have been optimised for use with Cellabs Positive Control Slide and Mounting Fluid. Evans Blue dye contained in the anti-mouse Ig-FITC reagent is a possible carcinogen, therefore avoid contact with the skin. Patient specimens and the positive control slide should be handled as though potentially infectious. A Positive Control Slide must be run with each test run. Refer to Material Safety Data Sheet (MSDS) for further information. For professional use only. To be used with confirmatory testing guidelines or standard quality assurance protocols for clinical samples.

SPECIMEN COLLECTION AND PREPARATION OF SLIDES

THROAT SWABS: Insert swab in the posterior nasopharynx and collect specimen. The swab may be smeared onto a slide and fixed in acetone or methanol (acetone preferred) for 5 minutes and allow to air dry. Another option is to place the swab in 0.5mL of Chlamydia transport medium (see below for recipe). Cut the swab shaft and leave swab submerged. Vortex the swab for 10 seconds and place 25µL on slide and fix in methanol for 5 minutes.

Chlamydia Transport Medium: 68.46g sucrose, 2.09g K₂HPO₄, 1.09g KH₂PO₄, Streptomycin 50 µg/mL, Vancomycin 100 µg/mL, and Fungizone 5 µg/mL. Dissolve sucrose, K₂HPO₄ and KH₂PO₄ in distilled water to 1 litre and autoclave at 115°C for 15 minutes. Add STERILE antibiotic solutions when cool. Store at or below 2-8°C.

NASOPHARYNGEAL ASPIRATES (NPA) & BRONCHOALVEOLAR LAVAGES (BAL): Place 20µL of unwashed NPA specimen on a single well slide, dry and fix in methanol or acetone for 5 minutes. Spin BAL or washed NPA specimens in a microfuge at 12000g for 5 minutes and discard the supernatant. Resuspend the pellet in 100 - 500 µL PBS or transport medium and place 20µL on slide. Dry and fix PBS suspensions in methanol or acetone for 5 minutes, and transport medium suspensions in methanol for 5 minutes. Fix cultures grown in plastic trays or bottles in methanol for 5 minutes, as acetone will affect the plastic.

INDUCED SPUTUM SAMPLES: Vortex 2-3mL sputum with the same volume of 0.3% dithiothreitol, then incubate for 5-10 minutes at 37°C. Incubate for 30-60 minutes if specimen is viscous. Add an equal volume of PBS, vortex, and then centrifuge at 12000g for 10 minutes (or 2000g for 30 minutes). Discard supernatant and resuspend pellet in 500 µL PBS. Place 25µL on slide and fix in methanol or acetone for 5 minutes.

If the fixed specimens are not tested immediately store at 2-8°C overnight or freeze at -20°C for up to 2 months.

INSTRUCTIONS FOR USE

1. Add 25µL of the to the fixed specimen and , covering well area.
2. Incubate the slides at 37°C in a moist chamber for 30 minutes. Do not allow the slides to dry as this may cause non-specific binding.
3. Rinse gently in a bath of PBS for about one minute.
4. Drain slide and remove excess moisture around the well with a tissue.
5. Add 25µL of the to the specimen and .
6. Repeat Steps 2 - 4.
7. Add a drop of to the slide well. Place a coverslip on top of the drop and remove air bubbles.
8. Scan the entire specimen well using a fluorescence microscope under oil immersion at X600-X1000 magnification. Read immediately or store at 2-8°C in the dark for up to 24 hours.

READING AND INTERPRETATION OF RESULTS

The most common chlamydial forms in specimens are extracellular EB's. These appear as bright apple-green fluorescent pin-point, smooth-edged disc shaped bodies (approximately 300nm in diameter) and can be seen against a background of reddish-brown counterstained cells. Reticulate bodies may also be observed. These are 2-3 times larger than the EB's and they either fluoresce evenly or possess dark centres with a halo of fluorescence. If inclusions are present their appearance will be similar to those in the positive control slide, i.e., a mass of EB's and RB's within an enclosing membrane adjacent to the host cell nucleus. Any material which can be distinguished from chlamydial forms or fluoresce other than apple-green should be disregarded. The control slide should be used for comparison with the appearance and size of EB's found in the specimen. A positive diagnosis can be made when fixed stained specimens show at least four or more chlamydial EB's. A negative diagnosis can be reported when fixed stained smears are free of chlamydial organisms but cells are present. These cells may be intact or ruptured columnar cells. Irregularly shaped fluorescent material that differs in size from chlamydial bodies described above or fluoresces white, red or yellow should be considered non-specific staining.

WASTE DISPOSAL

Dispose of any unused components as biohazardous waste. Where the anti-mouse Ig-FITC reagent has been disposed of in the sink, ensure it is flushed with large quantities of water (as the sodium azide it contains may react with copper/lead plumbing systems). For more information, please refer to the MSDS.

SENSITIVITY, SPECIFICITY, & OTHER DATA ON THE CHLAMYDIA CEL PN

Refer to summary table at end of insert. All data on the Chlamydia Cel Pn can be obtained in the product information sheet. Please ask your local distributor or contact Cellabs.

INDEMNITY NOTICE

Modifications or changes made in the recommended procedure may affect the stated or implied claims. A positive or negative result does not preclude the presence of other underlying causative agents. Cellabs and its agents and distributors shall not be liable for damages under these circumstances.

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